

FORUMTEC

Dental Solutions & Beyond

Shaping Futures with
Precision, Creativity, and
Advanced Technologies



About Forumtec

Three decades of innovation in dental technology

Founded in 1987 by a team of visionary engineers, Forumtec has emerged as a leader in electronic dental devices and advanced equipment, dedicated to transforming dental procedures with cutting-edge technology. From the groundbreaking Bingo-1020 apex locator to the innovative Statera vi cordless endo motor, our products are designed to ensure better, faster, and simpler dental treatments. Specializing in Endodontic and Implantology devices, we ensure quicker and precise treatments, enhancing both professional convenience and patient comfort.

Our legacy and vision

With over 30 years of experience, Forumtec continues to lead with innovation, redefining dental equipment and procedures that empower dental professionals with technology that leads to more accurate, efficient, and less invasive treatments. As we look to the future, we remain committed to pushing the boundaries of what's possible in dental technology.

Innovative solutions and cutting-edge technology, that empower dental professionals to achieve unparalleled success in their practices.



**Statera vi
Endo Motor**



**Wirele-X
Apex Locator**



**Spotter
Implant Detector**



**RomiApex A15
Apex Locator**

STATERA vi


Transform your endodontic practice with Statera vi where precision meets innovation



A Cordless Endo Motor with Integrated Apex Locator

- Exceptional versatility, allowing the use of a wide variety of files with the Balanced Force Adaptive Technique.
- Reciprocation motion with a 5° interval and adjustment range from 5° to 400°.
- Speed settings ranging from 50 to 3200 RPM and torque levels reaching up to 6.0 Ncm.
- High-resolution color screen.
- 360° Rotating contra-angle.
- Cradle with embedded power bank for wireless charging.
- Unique triangle shape designed for ergonomic comfort, ensuring stable placement on the tray.
- Brushless motor for a longer lifespan and reduced maintenance.

Forumtec's proprietary Balanced Force Adaptive Technique integrates the benefits of balanced reciprocation and rotational movement to optimize root canal preparation.

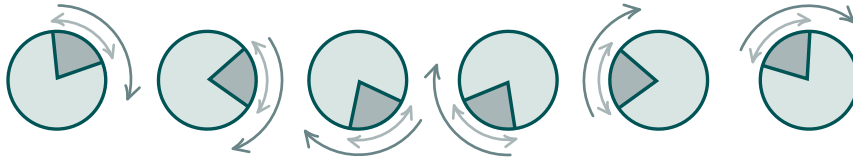

- The file performs a precise watch-winding motion, rotating back and forth with equal amplitudes, alternating between clockwise and counterclockwise rotations.
- This motion is repeated in several steps, increasing the amplitude with each cycle until the entire circumference of the canal is effectively covered.
- The system automatically adapts speed, torque, and the extent of the watch-winding rotation based on real-time canal anatomy and clinical conditions.



- Minimizes file separation, jamming, ledge formation, canal transportation, and perforation, even in narrow, curved canals.
- Lower apical extrusion than other reciprocating motions.
- Prevents self-threading and screw-in effects.
- Optimal control of axial and lateral movement.

wirelex

Upgrade your root canal treatments with the Wirele-X Apex Locator



Equipped with advanced Bluetooth connectivity for a seamless experience, its unique design and user-friendly interface deliver accurate measurements and smooth operation, making it an indispensable tool for modern dental practices

■ Innovative design

A compact apex locator that transmits measurements to a 7" display unit via Bluetooth, eliminating the complexity of handling long cables.

■ Enhanced user interface

A range of preferred working modes accessible through diverse graphic interfaces with customizable screen color choices.

■ Simplified measurement procedure

Demonstrates scale and numerical values along the root canal and beyond the apical foramen accompanied by audible alerts and numerical values.

■ Accurate and Reliable

High accuracy based on Forumtec's renowned patented measuring algorithm.

”

What do experts say about Wirele-X

“The accuracy and reliability of Wirele-X is simply outstanding. Combined with the convenience of having a wireless sensor, I am more confident in my measurements than ever before”

Nestor Cohenca, DDS | Prime Endodontics , US

“We have been using the Wirele-X Apex Locator in our Endodontic residency program at Albert Einstein Medical Center. We have found it to be very easy to use and extremely accurate. The wireless feature is most advantageous to have during endodontic treatment”

Fred Barnett, DMD | Einstein Medical Center, Philadelphia

“The Wirele-X is very accurate in the measurements in all clinical situations even in wet canals. The interface is very clear and gives a perfect view of the progression of the file in to the canal. It is very modern, easy to use and very reliable in every day practice”

Nicolas Gardon , DDS | French Society of Endodontics



spotter

Detect buried implants
accurately with Spotter
Implant detector



A minimal invasive, accurate, and time-saving solution for uncovering implants in a two-stage implantology procedure

■ **Accurate**

Detects buried implants beneath gum tissue or bone layer easily and accurately.

■ **Innovative**

Say goodbye to invasive procedures by eliminating the need to raise the flap.

■ **Cost-efficient**

Minimizes suturing requirements leading to a shorter and more comfortable recovery period eliminating complications.

■ **User-Friendly**

Compatible with all Titanium implant systems, offering flexibility and ensuring adaptability to different treatment scenarios.



Accurate detection of the center of buried implant



Ergonomic design for a more comfortable grip



Led indicators show the sensor's progress as it approaches the implant

1. Approaching the implant
2. Reverse direction
3. Detection of implant location



Detachable sensor heads for easy sterilization



The sensor holder can be adjusted to facilitate detection of implant in both upper and lower jaw



[Scan to watch a demo](#)

[Click to view a demo](#)

RomiApex A15

Discover the efficiency
of the RomiApex A15



Compact, powerful apex locator designed to maximize space and deliver premium features at an affordable price

■ Compact design

Lightweight and ergonomic design, optimizing space during treatment

■ Connection verification test

Includes cable connection verification test to ensure safe and reliable operation.

■ Automatic shutdown

The device automatically shuts down after 5 minutes, ensuring energy efficiency and safety.

■ Visual display feedback

Displays numerical values ranging from 2.0 to FO (Foramen Over) for clear measurement feedback.

■ Audio feedback

Audio feedback activates with file progression, offering four volume levels: mute, low, normal, and high.

Why choose Forumtec?



Legacy of Success

Over three decades of experience dedicated to continuous improvement and technological breakthroughs.



Innovative Solutions

Pioneering products that set new standards in dental care.



Trusted by Professionals

Holding a 35% market share in the Apex Locators category, distributed by top international brands.



Strong Client Partnerships


Commitment to building strong client partnerships, ensuring collaborative efforts and tailored solutions to meet their unique needs and goals.

FORUMTEC

Dental Solutions & Beyond

**Our commitment to innovation and precision allows us
to continuously redefine dental devices**

Visit our website

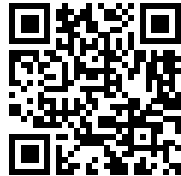

Click to visit our website

info@forumtec.net ■ www.forumtec.net
